

M498 Histosette® I


Tissue Processing/Embedding Cassettes

Made of acetal polymer

These cassettes are identical to Series M490 but **the anterior writing area has a 45° instead of a 30° angle**. The more acute angle makes these cassettes more suitable to be used with some models of cassette labelling instruments.

These disposable plastic cassettes hold tissue specimens very efficiently during the embedding process, as well as in a storage file. Moulded from a special high density polymer, these cassettes keep specimens safely submerged in liquid and are totally resistant to the chemical action of histological solvents. The efficient flow-through slots maximize fluid exchange and ensure proper drainage. The one-piece integral lid eliminates the need for separate steel lids. Just snap apart and conveniently lock the lid into the base of the cassette. They can be opened or closed as often as necessary and they always relock securely without danger of specimen loss. Large labelling areas are located on three sides of the cassette for your convenience.

Each case contains 3 dispenser boxes of 500 cassettes.

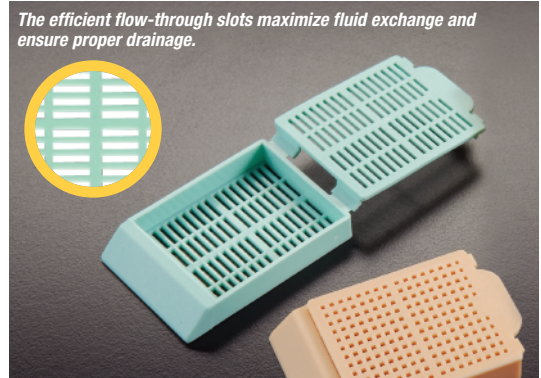

M499 Histosette® I

Biopsy Processing/Embedding Cassettes

Made of acetal polymer


This model of cassette is similar to Series M498 and is specially designed to **hold biopsy specimens during the embedding process**.

Each case contains 3 dispenser boxes of 500 cassettes.



	tissue code	biopsy code	colour	units
Tissue/Biopsy Processing/Embedding Cassettes	M498-2	M499-2	White	1500
Tissue/Biopsy Processing/Embedding Cassettes	M498-3	M499-3	Pink	1500
Tissue/Biopsy Processing/Embedding Cassettes	M498-4	M499-4	Green	1500
Tissue/Biopsy Processing/Embedding Cassettes	M498-5	M499-5	Yellow	1500
Tissue/Biopsy Processing/Embedding Cassettes	M498-6	M499-6	Blue	1500
Tissue/Biopsy Processing/Embedding Cassettes	M498-7	M499-7	Peach	1500
Tissue/Biopsy Processing/Embedding Cassettes	M498-8	M499-8	Tan	1500
Tissue/Biopsy Processing/Embedding Cassettes	M498-9	M499-9	Grey	1500
Tissue/Biopsy Processing/Embedding Cassettes	M498-10	M499-10	Lilac	1500
Tissue/Biopsy Processing/Embedding Cassettes	M498-11	M499-11	Orange	1500
Tissue/Biopsy Processing/Embedding Cassettes	M498-12	M499-12	Aqua	1500
Tissue/Biopsy Processing/Embedding Cassettes	M498-16	M499-16	Fluo Pink	1500
Tissue/Biopsy Processing/Embedding Cassettes	M498-17	M499-17	Fluo Green	1500
Tissue/Biopsy Processing/Embedding Cassettes	M498-18	M499-18	Fluo Yellow	1500

Anatomy of a Histosette® I

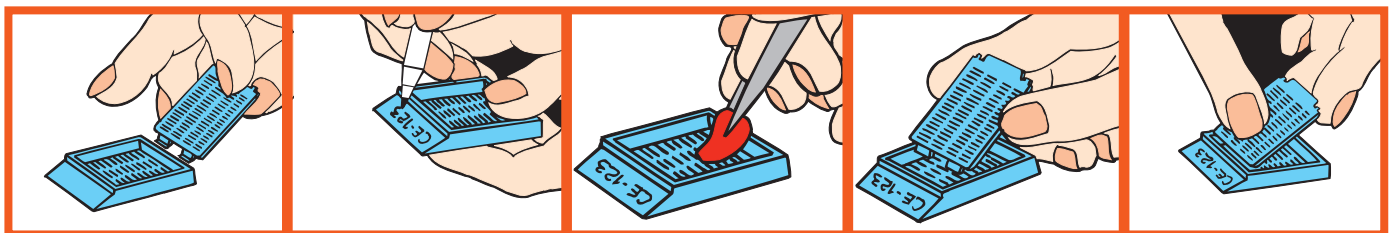

Separate cover by folding it forward and backward.

Identify in front or sides, using M495-10 Marker Pen.

Insert sample.

Slide cover in place and lock it.

After processing, remove cover by pulling on back tab.